

Med-Lux



How it works

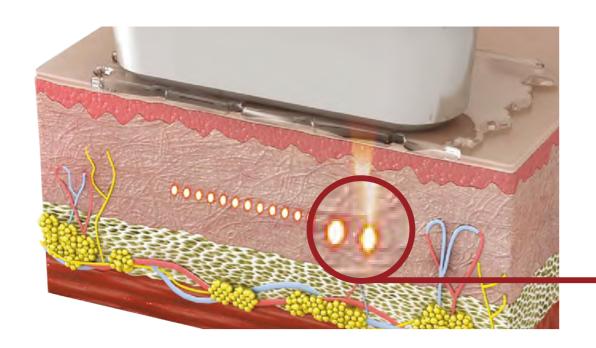

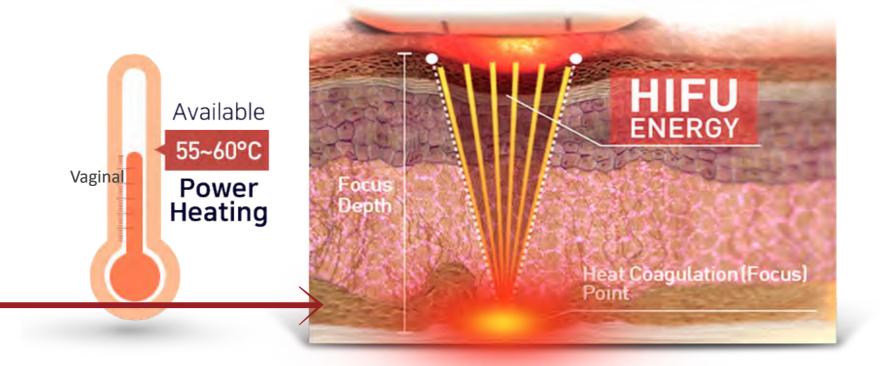

HIFU-powered transducers deliver concentrated waves of ultrasound into multiple layers of skin and body tissues without making direct contact with surrounding areas.

Transferring stabilized temperatures of 65~75°C into the targeted zones, the SMASFOCUS V aims to induce a process of neocollagenesis in the face while tightening body tissues for contouring purposes


FOR LIFTING & TIGHTENING



1 The dermis layer of the skin naturally synthesizes less collagen over the course of time.



3 The dermis layer of the skin naturally synthesizes less collagen over the course of time.

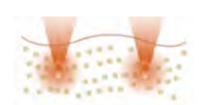

2 Focused ultrasound waves deliver energy directly within the skin layer causing a shrinking effect.



Through the regeneration of collagen, wrinkles are eliminated in the process of both firming and tightening the skin.


FOR SLIMMING & CONTOURING



1 The dermis layer of the skin naturally synthesizes less collagen over the course of time



3 As a natural healing response of the body the cells regenerate joining together to fill the gap.

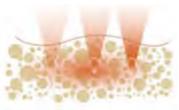



2 The SMASFOCUS V energy destroys the fat cells that have been targeted.



4 In result, the layer of fat naturally exits the body delivering a smooth, contoured body line.


FOR VAGINAL TIGHTENING

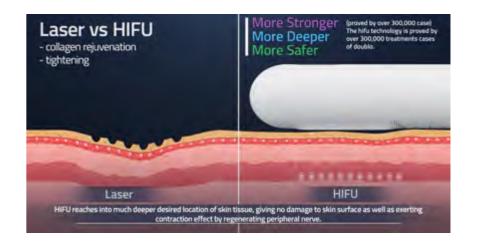

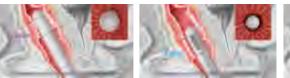



Product Introduction

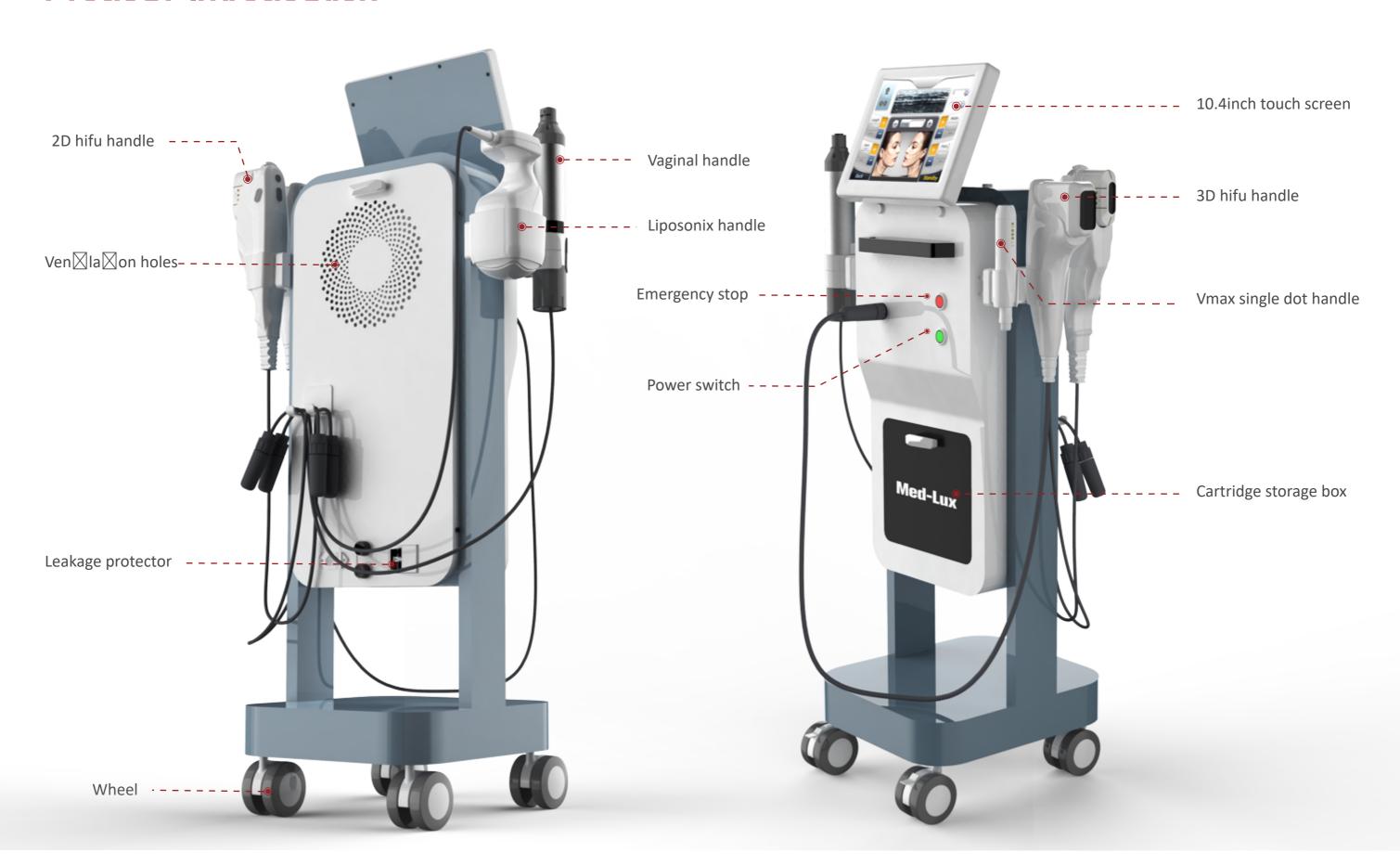



Vmax single dot handle



1.5mm, 3.0mm, 4.5mm





2D hifu handle



1.5mm, 3.0mm, 4.5mm

Check, Forehead, Double chin, Neck wrinkles, V-line forming



3D hifu handle




1.5mm, 3.0mm, 4.5mm



6.0mm, 8.0mm, 10.0mm, 13.0mm, 16.0mm

Check, Forehead, Double chin, Neck wrinkles, V-line forming

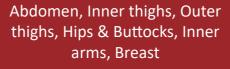



Liposonix handle



1.3 cm

8.0mm, 13mm

Body slimming and fat reduction



Vaginal handle



3.0mm, 4.5mm

Vaginal tightening











Easy to upgrade


Each machine is pre-loaded with a software system with all handle functions before shipment. When users need to upgrade the machine functions, they only need to purchase the handle with that function, plug the handle into the machine, and the software system can automatically identify the new handle to realize machine function upgrade.



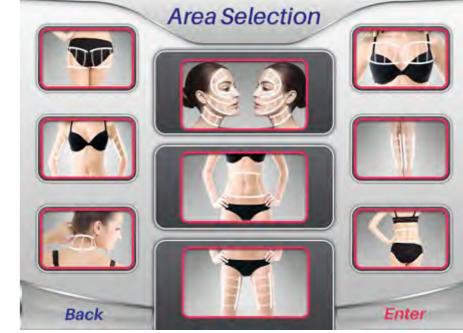

2D HIFU 3D HIFU

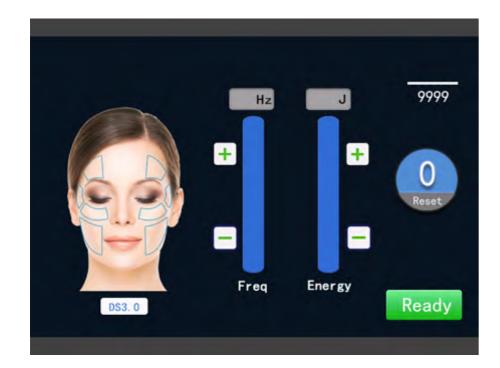

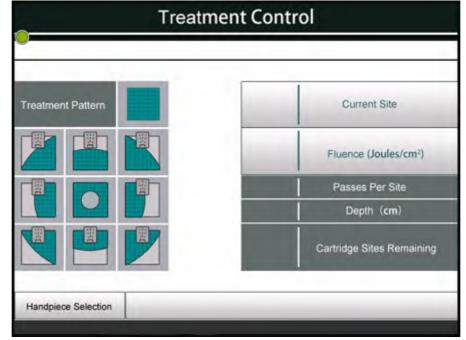



Vmax Liposonix Vaginal


Summary of Advantages

- Procedure Time 30 45 Minutes
- Full recovery Immediate
- Number of Treatments

 1 every 12 18 months

- Back To Work Immediate
- Software Perloaded
 All func⊠ons' so⊠ware
 perloaded already
- Risks & Complications

 Possible mild redness,
 swellling and tenderness

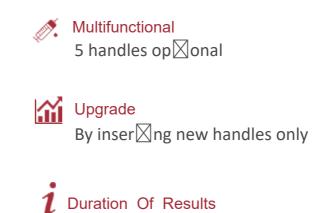

12 - 18 months



Before and after



Before After Before After Before After

Specifications

Power		1200 W
Frequency	1	4 ~ 7MHz
Net Weigh	t	25kg
Dimension		400*350*1060mm
Voltage AC		100-240v 50/60Hz

